



Version 2.0706

**EpiQuik™ Global Histone H3-K9
Methylation Assay Kit**

Catalog No. P-3018

User Guide*

***Always use the most updated User
Guide included in your current order.**

EpiGenetek Group Inc.
760 Parkside Avenue
Brooklyn, NY 11226
Tel: 1 877-374-4368
Fax: 1 718-484-3956
info@epigenetek.com

www . EPIGENETEK . com

FOR RESEARCH USE ONLY. NOT FOR USE IN DIAGNOSTIC PROCEDURES.

TABLE OF CONTENTS

Introduction	3
Principle and Procedure	4
Product Use Information	5
Kit Contents	6
Shipping and Storage	6
Materials Required But Not Supplied	6
Protocol	7
Troubleshooting	10
Ordering Information	11

Epigentek Group Inc.
760 Parkside Avenue
Brooklyn, NY 11226
Tel: 1 877-374-4368
Fax: 1 718-484-3956
info@epigentek.com


INTRODUCTION

Epigenetic activation or inactivation of genes play a critical role in many important human diseases, especially in cancer. A major mechanism for epigenetic inactivation of the genes is methylation of CpG islands in genome DNA caused by DNA methyltransferases. Histone methyltransferases (HMTs) control or regulate DNA methylation through chromatin-dependent transcription repression or activation. HMTs transfer 1-3 methyl groups from S-adenosyl-L-methionine to the lysine and arginine residues of histone proteins. ESET, G9a, SUV39-h1, SUV39-h2, SETDB1, Dim-5, and Eu-HMTase, are histone methyltransferases that catalyze methylation of histone H3 at lysine 9 (H3-K9) in mammalian cells. H3-K9 methylation mediates heterochromatin formation by forming a binding site for HP1 and also participates in silencing gene expression at euchromatic sites. Increased global H3-K9 methylation is also found to be involved in some pathological processes such as cancer progress. There are only a couple of methods such as western blot used for measuring histone H3-K9 methylation. These methods available so far are time consuming and labor intensive, or have low throughput. The *EpiQuik*[™] Global H3-K9 Methylation Assay Kit addresses these problems by using a unique procedure to measure global methylation of histone H3-K9. The kit has the following features:

- Quick and efficient procedure, which can be finished within 5 hours.
- Innovative colorimetric assay without the need for radioactivity, electrophoresis, or chromatography.
- Strip microplate format makes the assay flexible: manual or high throughput.

PRINCIPLE AND PROCEDURE

The *EpiQuik*™ Global H3-K9 Methylation Assay Kit is designed for measuring global histone H3-K9 methylation. In an assay with this kit, the histone proteins are stably spotted on the strip wells. The methylated histone H3-K9 can be recognized with a high-affinity antibody. The ratio or amount of methylated H3-K9 can be quantified through HRP conjugated secondary antibody-color development system and is proportional to the intensity of color development.

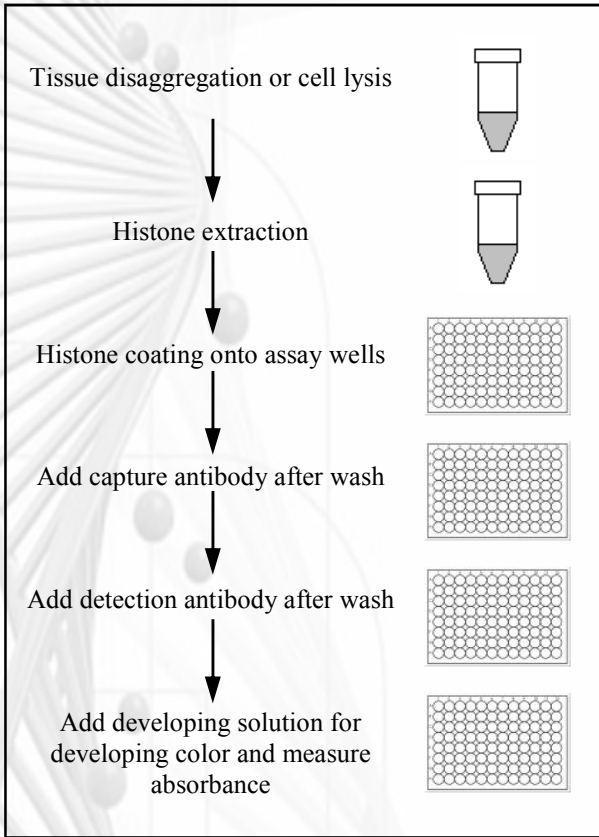

Schematic Procedure for Using the *EpiQuik*™ Global Histone H3-K9 Methylation Assay Kit

PRODUCT USE INFORMATION

The *EpiQuik*[™] Global Histone H3-K9 Methylation Assay Kit is suitable for specifically measuring global histone H3-K9 methylation using a variety of mammalian cells including fresh and frozen tissues, and cultured adherent and suspension cells.

The *EpiQuik*[™] Global Histone H3-K9 Methylation Assay Kit for research use only and is not intended for diagnostic or therapeutic application.

Suitable lab coat, disposable gloves, and eye protection is required when working with the kit.

Epigentek guarantees the performance of all products in the manner described in our product instructions.

Epigentek reserves the right to change or modify any product to enhance its performance and design.

EpiQuik[™] is a trademark of Epigentek, Inc.

The *EpiQuik*[™] Global Histone H3-K9 Methylation Assay Kit and methods of use are covered by a pending US patent.

KIT CONTENTS

Components	48 assays P-3018-48	96 assays P-3018-96
GC1 (10X Lysis Buffer)	5 ml	10 ml
GC2 (Extraction Buffer)	8 ml	16 ml
GC3 (10X Wash Buffer)	14 ml	28 ml
GC4 (Histone Buffer)	0.5 ml	1 ml
GC5 (Blocking Buffer)	10 ml	20 ml
GC6 (Antibody Buffer)	6 ml	12 ml
GC7 (Capture Antibody, 100 µg/ml)*	25 µl	50 µl
GC8 (Detection Antibody, 400 µg/ml)*	10 µl	20 µl
GC9 (Developing Solution)	5 ml	10 ml
GC10 (Stop Solution)	3 ml	6 ml
Methylated H3-K9 Control (60 µg/ml)*	10 µl	20 µl
8-Well Assay Strips (with frame)	6	12
User Guide	1	1

* For maximum recovery of the products, centrifuge the original vial prior to opening the cap.

SHIPPING AND STORAGE

The kit is shipped in two parts: one part at ambient room temperature, and the second part on frozen ice packs at 4°C.

Upon receipt: (1) Store **GC8** and the **Methylated H3-K9 Control** at -20°C; (2) Store **GC3**, **GC5**, **GC7**, **GC9**, and **8-Well Assay Strips** at 4°C away from light; (3) Store **all other components** at room temperature. The kit is stable for up to 6 months from the shipment date, when stored properly.

Note: Check if wash buffer, **GC3**, contains salt precipitates before using. If so, warm (at room temperature or 37°C) and shake the buffer until the salts are re-dissolved.

MATERIALS REQUIRED BUT NOT SUPPLIED

Orbital shaker	Dounce homogenizer
Pipettes and pipette tips	100 % TCA solution
Microplate reader	Glycerol
1.5 ml microcentrifuge tubes	Acetone
60 or 100 mm plate	5% HCl

PROTOCOL

Nucleic Extraction Preparation

For Tissue Samples:

1. Place the tissue sample into a 60 or 100 mm plate. Remove unwanted tissue, such as fat and necrotic material from the sample. Weigh the sample and cut the sample into small pieces (1-2 mm³) with a scalpel or scissors.
2. Transfer tissue pieces to a Dounce homogenizer. Dilute **GC1** with distilled water (pH 7.2 to 7.5) at a 1:10 ratio (e.g., 1 ml of **GC1** + 9 ml of distilled water). Add 1 ml of the **diluted GC1** per every 200 mg of tissue and disaggregate tissue pieces by 10-30 strokes.
3. Transfer homogenized mixture to a 15 ml conical tube and centrifuge at 3000 rpm for 5 minutes at 4°C. If total mixture volume is less than 2 ml, transfer mixture to a 2 ml vial and centrifuge at 10,000 rpm for 1 minute at 4°C. Remove supernatant.

For Adherent Cells:

1. Cells (treated or untreated) are grown to 70-80% confluency, then trypsinized and collected into a 15 ml conical tube. Count cells in a hemacytometer.
2. Centrifuge the cells at 1000 rpm for 5 minutes and discard the supernatant. Wash cells with 10 ml of PBS once by centrifugation at 1000 rpm for 5 minutes. Discard the supernatant.
3. Add the **diluted GC1** to re-suspend cell pellet (200 μ l/1 x 10⁶ cells). Transfer cell suspension to a 1.5 ml vial and incubate on ice for 5 minutes and vortex occasionally.
4. Pellet cell debris by centrifuging at 12,000 rpm for 30 seconds.

For Suspension Cells:

1. Collect cells (treated or untreated) into a 15 ml conical tube. (1-2 x 10⁶ cells are required for each reaction). Count cells in a hemacytometer.
2. Centrifuge the cells at 1000 rpm for 5 minutes and discard the supernatant. Wash cells with 10 ml of PBS once by centrifugation at 1000 rpm for 5 minutes. Discard the supernatant.
3. Add the **diluted GC1** to re-suspend cell pellet (100 μ l/1 x 10⁶ cells). Transfer cell suspension to a 1.5 ml vial and incubate on ice for 5 minutes and vortex occasionally.
4. Pellet cell debris by centrifuging at 12,000 rpm for 30 seconds. Discard the supernatant.

(Continued on Next Page)

Histone Extraction

1. Add *Glycerol* to **GC2** at a 1:10 ratio (ex: add 100 μl of *Glycerol* to 900 μl of **GC2**) to prepare the **GC2/Glycerol Solution**. Add the **diluted GC1** to cell debris (10 μl /1 $\times 10^6$ cells or 40 mg of tissue), followed by adding 3 volumes of **GC2/Glycerol Solution**. Mix by vortexing and incubate on ice for 5 minutes.
2. Pellet nucleic debris by centrifuging at 12,000 rpm for 5 minutes at 4°C. Transfer the supernatant to a 1.5 ml vial.
3. Add **100% TCA Solution** to the supernatant at a 1:4 ratio (ex: add 100 μl of TCA to 300 μl of supernatant; final concentration of TCA should be 25%). Incubate on ice for 30 minutes.
4. Collect the precipitate by centrifuging at 12,000 rpm for 10 minutes at 4°C.
5. Remove supernatant and add 1 ml of *acetone* containing 0.1% HCl to precipitate. Mix and incubate on ice for 1 minute.
6. Collect the pellet by centrifuging at 12,000 rpm for 2 minutes at 4°C. Wash the pellet with 1 ml of undiluted *acetone*. Allow 1 minute on ice for wash. Discard the supernatant.
7. Collect the pellet by centrifuging at 12,000 rpm for 2 minutes at 4°C. Remove the supernatant as much as possible and air dry the pellet for 5 minutes.
8. Add distilled water to dissolve pellet (10 μl of water per amount of pellet extracted from 1 $\times 10^6$ cells or 40 mg of tissue) and measure histone protein concentration. The histone extract can be used immediately or stored at -20°C.

Histone H3-K9 Methylation Detection

1. Determine the number of strip wells required. Leave these strips in the plate frame (remaining unused strips can be placed back in the bag. Seal the bag tightly and store at 4°C). Dilute **GC3** with distilled water (pH 7.2-7.5) at a 1:10 ratio.
2. Adjust protein concentration to 200 ng/ μl or 400 ng/ μl with **GC4** and add 5 μl (1-2 μg) of the protein solution into the central area of each well. Spread out the solution over the bottom of the strip well by pipetting the solution up and down several times. Incubate the wells at 37°C (with no humidity) for 90 minutes to evaporate the solution and dry the wells. For the blank, add 5 μl of **GC4** to the wells. For the positive control, dilute the **Methylated H3-K9 Control** to 2-30 ng/ μl with **GC4**, and add 5 μl (10-150 ng) of the **diluted Methylated H3-K9 Control** solution to the wells.

3. Add 150 μl of **GC5** to the dried wells and incubate at 37°C for 30-45 minutes.
4. Aspirate and wash the wells with 150 μl of the **diluted GC3** three times.
5. Dilute **GC7** (at a 1:100 ratio) to 1 $\mu\text{g}/\text{ml}$ with **GC6**. Add 50 μl of the **diluted GC7** to the wells and incubate at room temperature for 60 minutes on an orbital shaker (50-100 rpm).
6. Aspirate and wash the wells with 150 μl of the **diluted GC3** four times.
7. Dilute **GC8** (at a 1:1000 ratio) to 0.4 $\mu\text{g}/\text{ml}$ with **GC6**. Add 50 μl of the **diluted GC8** to the wells and incubate at room temperature for 30 minutes.
8. Aspirate and wash the wells with 150 μl of the **diluted GC3** five times.
9. Add 100 μl of **GC9** to the wells and incubate at room temperature for 2-10 minutes away from light. Monitor the color development in the sample and standard wells (blue).
10. Add 50 μl of **GC10** to the wells and read absorbance on a microplate reader at 450 nm.
11. Calculate % H3-K9 methylation:

$$\text{Methylation \%} = \frac{\text{OD (sample - blank)}}{\text{OD (untreated control - blank)}} \times 100\%$$

For accurate calculation, plot OD value versus amount of **Methylated H3-K9 Control** and determine the slope as $\Delta \text{OD}/\text{ng}$.

Calculate the amount of methylated H3-K9 using the following formula:

$$\text{Amount (ng/mg protein)} = \frac{\text{OD (sample - blank)}}{\text{slope}} \times 1000$$

TROUBLESHOOTING

No Signal for Both the Positive Control and the Samples

Reagents are added incorrectly. Check if reagents are added in proper order and if any steps of the procedure may have been omitted by mistake.

The well is not completely dried. Ensure the well is incubated with no humidity and dry before adding the Blocking Buffer.

The well is incorrectly washed before protein coating. Ensure the well is not washed before adding the positive control or protein extracts.

Incubation time and temperature is incorrect. Ensure the incubation time and temperature described in the protocol are followed correctly.

No Signal or Very Weak Signal for Only the Positive Control

The positive control is added insufficiently. Ensure a sufficient amount of control is added to the well.

No Signal for Only the Sample

The protein sample is not properly extracted. Follow the protocol instructions for the histone protein extraction.

The protein amount is added into well insufficiently. Ensure extract contains a sufficient amount of protein.

Protein extracts are incorrectly stored or have been stored for a long time. Ensure the protein extracts are stored at -20°C for no more than 8 weeks.

High Background Present for the Blank

The well is not washed sufficiently. Check if wash at each step is performed according to the protocol.

Contaminated by the positive control. Ensure the well is not contaminated from adding the control protein or from using control protein contaminated tips.

Overdevelopment. Decrease development time in Step 9 of "Histone H3-K9 Methylation Detection."

ORDERING INFORMATION

Products	Size	Cat. No.
<i>EpiQuik</i> [™] Global Histone H3-K9 Methylation Assay Kit	48 assays	P-3018-48
	96 assays	P-3018-96

Available Related Products

EpiQuik[™] Histone Methyltransferase Activity/Inhibition Assay Kit (H3-K4) P-3002

EpiQuik[™] Histone Methyltransferase Activity/Inhibition Assay Kit (H3-K9) P-3003

EpiQuik[™] *In Situ* Histone H3-K4 Methylation Assay Kit P-3015

EpiQuik[™] *In Situ* Histone H3-K9 Methylation Assay Kit P-3016

EpiQuik[™] Global Histone H3-K4 Methylation Assay Kit P-3017

Need more components? You can also order parts separately by calling 1-877-374-4368 or e-mailing sales@epigentek.com.

