



pBABEpuro-hTERT Retroviral Vector

CATALOG NUMBER: RTV-006

STORAGE: -80°C

QUANTITY AND CONCENTRATION: 100 µL of bacterial glycerol stock

Background

Retroviruses are efficient tools for delivering heritable genes into the genome of dividing cells. Cell Biolabs' retrovirus vector is based on the pBABE vector system, which is derived from Moloney murine leukemia virus (MMLV). The vector provides the viral package signal, transcription and processing elements, and a target gene. The viral *env* gene, produced by the package cell line, encodes the envelop protein, which determines the viral infectivity range. Transfection into a package cell line produces high-titer, replication-incompetent viruses. In addition to transfer and expression of exogenous genes in mammalian cells, recently, retroviruses are used to express silencing RNAs (siRNA) to decrease the expression of target genes both *in vitro* and *in vivo*.

The vector contains the bacterial origin of replication, ampicillin-resistance gene, and Puromycin-resistance gene for the growth of infected mammalian cells to select stable cell lines (Figure 1).

hTERT, or "human Telomerase Reverse Transcriptase," is a ribonucleoprotein that maintains telomere ends by addition of the telomere repeat sequence TTAGGG. These TTAGGG repeats (with their various protein binding partners) are called telomeres. Telomerase expression plays a role in cellular senescence, as it is normally repressed in postnatal somatic cells resulting in progressive shortening of telomeres. Deregulation of telomerase expression in somatic cells may be involved in oncogenesis. A hTERT gene with an HA tag at its C-terminal is cloned into the retroviral vector pBABEpuro between EcoRI and Sall sites.

Safety Consideration

Remember that you will be working with samples containing infectious virus. Follow the recommended NIH guidelines for all materials containing BSL-2 organisms. Always wear gloves, use filtered tips and work under a biosafety hood.

References

1. Morgenstern, J. P. and H Land. (1990) *Nuc. Acid Res.* 18, 3587-3596.
2. Coffin, J. M. and H. E. Varmus, *Retroviruses*, Cold Spring Harbor Press, NY.
3. Schuck S, Manninen A, Honsho M, Fullekrug J and Simons K. (2004) *Proc Natl Acad Sci U S A.* 101, 4912-4917.

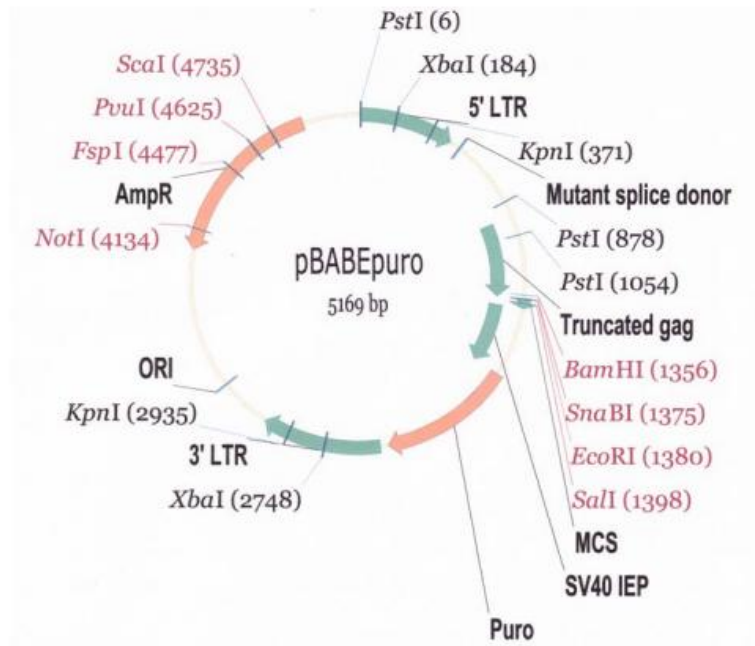


Figure 1. pBABE-Puro Retroviral Vector Map

Warranty

These products are warranted to perform as described in their labeling and in Cell Biolabs literature when used in accordance with their instructions. THERE ARE NO WARRANTIES THAT EXTEND BEYOND THIS EXPRESSED WARRANTY AND CELL BIOLABS DISCLAIMS ANY IMPLIED WARRANTY OF MERCHANTABILITY OR WARRANTY OF FITNESS FOR PARTICULAR PURPOSE. CELL BIOLABS's sole obligation and purchaser's exclusive remedy for breach of this warranty shall be, at the option of CELL BIOLABS, to repair or replace the products. In no event shall CELL BIOLABS be liable for any proximate, incidental or consequential damages in connection with the products.

This product is for RESEARCH USE ONLY; not for use in diagnostic procedures.

GEL